

Dragon Professional Anywhere 2025.3

Audio Routing Guide

Trademarks

Nuance®, the Nuance logo, Dragon, the Dragon logo, PowerScribe™ and PowerMic™ are trademarks or registered trademarks of Nuance Communications, Inc. and/or its subsidiaries in the United States and/or other countries.

Other names and trademarks referenced herein are trademarks or registered trademarks of their respective owners.

Copyright notice

This publication is copyrighted, and all rights are reserved by Nuance Communications, Inc. No part of this publication may be reproduced, transmitted, transcribed, stored in a retrieval system, or translated into any language or computer language, in any form or by any means, electronic, mechanical, magnetic, optical, chemical, manual, or otherwise, without the prior written permission of Nuance Communications, Inc., 1 Wayside Rd., Burlington, MA 01803.

Copyright © 2002-2025 Nuance Communications, Inc. All rights reserved.

Disclaimer

All content within this guide is intended for general information only. Nuance, its affiliates, agents and employees, and all persons acting on its or their behalf ("Nuance"), disclaims to the maximum extent permitted by law, any and all liability that may be based on this information or errors or omissions in this guide. Nothing stated or posted in this guide or available through any services is intended to be, and must not be taken to be, the practice of medicine or the provision of medical care.

Reasonable efforts were used to present an accurate guide. Given the constant advances in technology, medical knowledge, and legal and regulatory frameworks in which healthcare services operate, Nuance reserves the right to change, amend, update or replace the information contained in this guide, in part or entirely, for any reason or no reason. Nuance will endeavor, without implied guarantee, to keep the content and information contained in this guide reasonably current and accessible, but Nuance will not be liable for environmental and operational changes, errors or omissions, or factors it does not control, including without limitation, failure of electronic and mechanical equipment, communication over the internet, computer viruses, unauthorized access or misuse. Nuance does not guarantee, represent, or warrant that your use of the services will be uninterrupted, error-free, accurate, reliable or will meet your requirements. Nuance does not guarantee that you will be able to access or use the platform (either directly or through third-party networks) at the time or location of your choosing. Mention of non-Nuance products or services is for information purposes only and constitutes neither an endorsement nor a recommendation.

To the maximum extent permitted by law, the information in this guide is provided "AS IS," "AS AVAILABLE," "WITH ALL FAULTS," and without warranty of any kind, either express or implied, including but not limited to, the implied warranties of merchantability, fitness for a particular

purpose, and non-infringement. Nuance shall not under any circumstances be liable to any person for any special, incidental, indirect or consequential damages, including, without limitation, damages resulting from use of or reliance on the information presented, loss of profits or revenues or costs of replacement goods, even if informed in advance of the possibility of such damages. Your use of, or inability to use, or activity in connection with the services and information in this guide is at your sole risk. Nuance expressly disclaims any responsibility for any damage, injury, harm, cost, expense, or liability arising out of or related to your use or misuse of this guide.

Nuance Communications
One Wayside Road
Burlington, MA 01803
www.nuance.com/healthcare

Table of contents

Overview	5
System configurations	6
Single-hop virtualization	6
Double-hop virtualization	7
Available solutions	9
Nuance virtual extensions	9
Vendor native audio channels	10
USB redirection	11
Recommendations	12
Zero clients	12
Single-hop configuration (thin or thick clients)	12
Double-hop configuration (thin or thick clients)	12
Recommendations overview	13
References	14

Overview

For accurate speech recognition, Dragon Professional Anywhere must capture clear audio from the user's device. To achieve this in a virtualized environment, it requires the following considerations:

- The client workstation, where the user is dictating, is not running a locally-installed instance of Dragon Professional Anywhere or the target app. The apps are running on a virtualization server or virtual desktop. The workstation displays a bitmap representation of the apps or virtual desktop via a receiver application.
- Dragon Professional Anywhere must be installed on the same server or virtual desktop as the target app to access the app's text fields.
- Button controls and high-quality audio must be routed from the workstation where the PowerMic is connected to the virtualization layer where Dragon Professional Anywhere and the target app are running.

This guide summarizes the types of virtualized environments at customer sites and the available Nuance solutions. It helps Nuance customer-facing personnel to recommend the best solution for each customer's configuration.

Note: Many sites install Dragon Professional Anywhere locally in addition to a virtual configuration; local installations aren't the primary focus of this document.

System configurations

Dragon Professional Anywhere and the target app can be hosted on a virtualization server (for example, Citrix XenApp) or included in a virtual desktop image running on a virtualization server (for example, Citrix XenDesktop).


The workstation can be a thick client running a full Windows version, a thin client running Linux or Windows Embedded or a zero client with no operating system.

A site's configuration might include a virtualization vendor for which Nuance provides custom channels (Citrix, Omnisia/VMware, Microsoft RDS/AVD); this affects the choice of solution for the site.

A customer might use more than one of these configurations. Therefore, more than one Nuance solution might be needed to fulfil a site's requirements.


Single-hop virtualization

Dragon Professional Anywhere and the target app are hosted on a virtualization server or included in a virtual desktop. The apps/virtual desktop are streamed to the user's workstation. Audio and microphone buttons must be routed from the workstation to the virtualization server/virtual desktop where Dragon Professional Anywhere and the target app are installed.

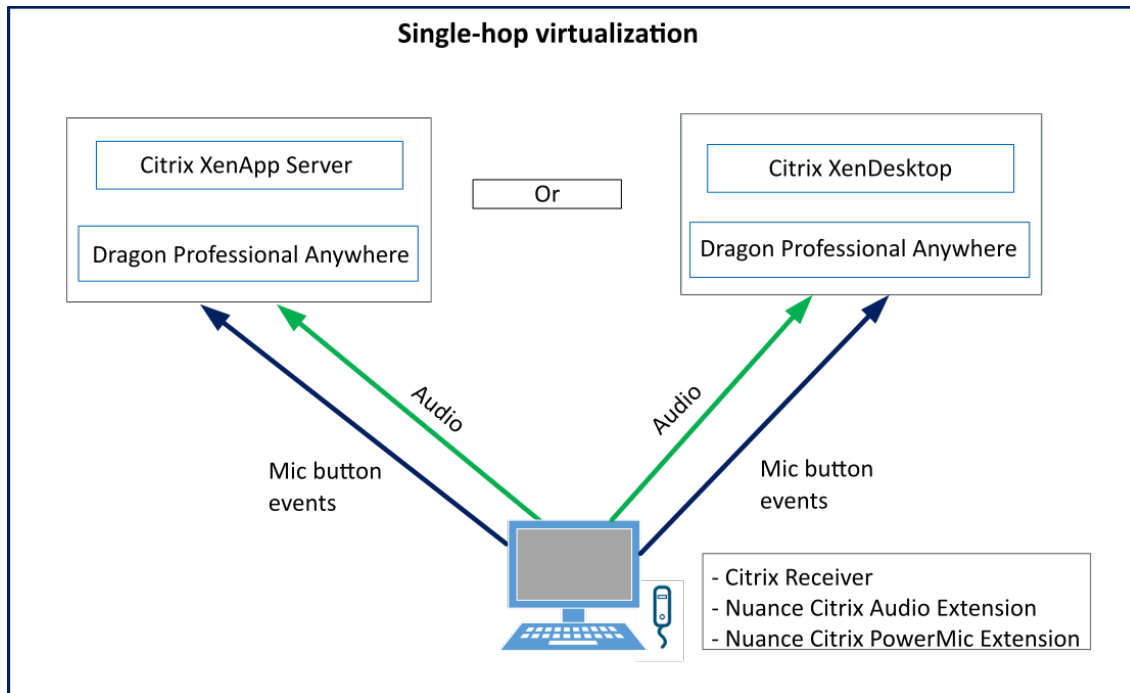

Double-hop virtualization

A double-hop configuration is most common when a site accesses an application hosted by a vendor and also uses a desktop virtualization solution such as XenDesktop. In this configuration, Dragon Professional Anywhere is hosted with the target app. Audio and microphone buttons must be routed through two virtualization layers or delivered via an alternative route that bypasses the intermediate virtualization layer.

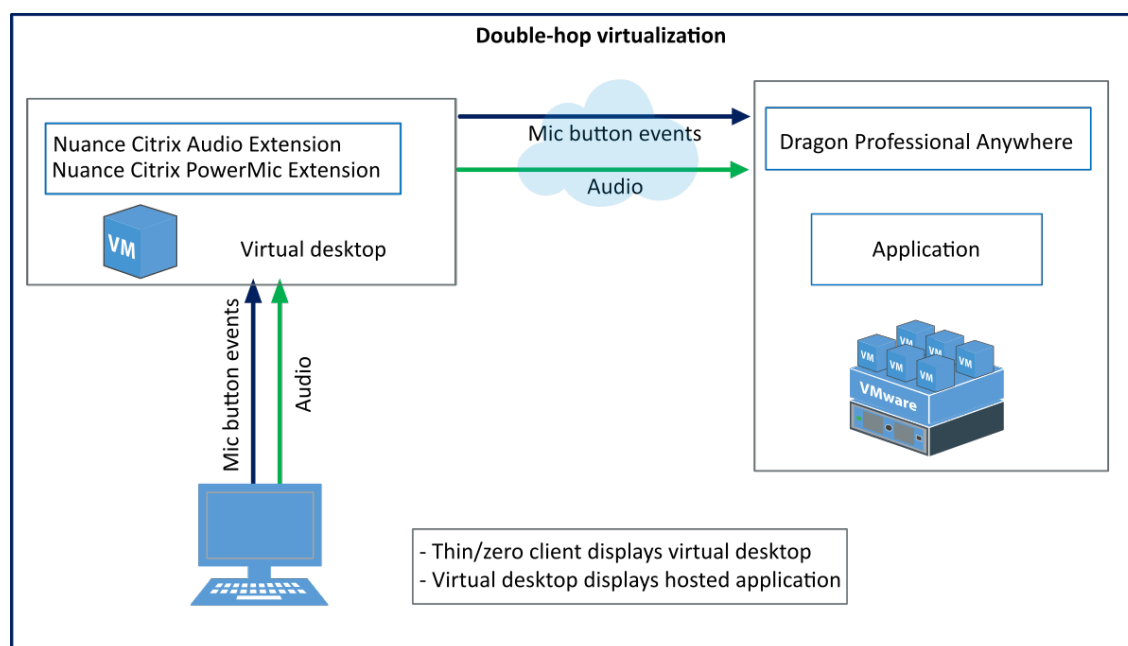

For more information on double-hop configuration for Citrix environments, see *Nuance - Double Hop Configuration for Citrix Environments* included in the Nuance virtual extensions package.

Available solutions

Nuance virtual extensions

The Nuance Citrix, Omnissa (formerly VMware) and RDS audio extensions and PowerMic extensions provide custom audio and microphone button channels for the following products:

- Citrix XenApp and XenDesktop/Citrix Virtual Apps and Desktops
- Remote Desktop Services (RDS)
- Azure Virtual Desktop (AVD)
- Omnissa Horizon
- Microsoft Windows and Windows Embedded operating systems
- Linux thin clients (for more information on thin clients supporting device splitting, contact your vendor)
- The following speech recognition apps:

Dragon Professional Anywhere

Apps based on Dragon Medical SpeechKit (.NET and COM editions)

Audio extensions

The Nuance Citrix/Omnissa/RDS audio extensions reduce the bandwidth required to transfer audio data from the client workstation to the server/virtual desktop where Dragon Professional Anywhere is installed. Native audio channels can require up to 1.4 Mbit/s bandwidth between the client workstation and the hosted app. The Nuance audio extensions require 27 kbit/s. This enables Nuance to offer compatibility with the most common virtualization solutions while at the same time delivering industry-leading speech recognition accuracy.

The Nuance audio extension is only installed on the client workstation (client end point).

PowerMic extensions

Microphone buttons, sliders and other controls must be routed to Dragon Professional Anywhere separately from audio. To enable this for the PowerMic in a Citrix XenApp/XenDesktop, Omnissa or Microsoft RDS/AVD system, the Nuance PowerMic extensions provide a custom channel for PowerMic button controls.

The Nuance PowerMic extension is only installed on the client workstation (client end point).

Azure Virtual Desktop

The Nuance RDS Client Audio Extension and Nuance PowerMic RDS Client Extension support the following clients:

- Windows Desktop client (MSRDC) for Azure Virtual Desktop.
- Windows App for access to Windows 365, Azure Virtual Desktop, and Microsoft Dev Box.
- Remote Desktop client (MSTSC) for Windows, Mac, iOS and Android.

These extensions don't support Remote Desktop and AVD clients for Web or clients installed from the Microsoft Store app.

Vendor native audio channels

Important: To use PowerMic buttons, you need to configure device splitting (available for Citrix and Omnisia/VMware environments, Windows client end points and some Linux-based clients). For more information on thin clients supporting device splitting, contact your vendor.

Citrix HDX

Citrix HDX audio channel (Citrix optimized virtual channel) redirects audio data to the virtual desktop at a significantly lower bandwidth than regular USB redirection.

Microsoft RDP

For more information about the bandwidth specification for Remote Desktop Protocol (RDP), see the [Microsoft documentation](#).

Omnissa RTAV

Omnissa (formerly VMware) Real Time Audio and Video (RTAV) redirects audio and video data to the virtual desktop at a significantly lower bandwidth than regular USB redirection.

Limitations

If you use native audio channels instead of the Nuance audio extensions, you might experience the following limitations:

- The bandwidth increases from 28 kbit/s to 150 kbit/s or higher.
- The microphone input volume can't be changed on the client end point; this might lead to audio gaps if the volume isn't configured correctly.
- Automatic gain control (AGC) and voice activity detection (VAD) aren't available.
- Controlling the client end point isn't possible (for example, disabling the screen saver during recording or disabling the buttons/stopping recording when the screen is locked).
- In standby mode, the audio data are streamed continuously from the client end point to the server or virtual desktop. When recording is on, audio data are also streamed even if the user isn't speaking.
- If Citrix HDX is used in combination with USB device splitting, each time a user roams to another client end point, they must manually select the HID device they want to redirect to be able to use the microphone buttons.

USB redirection

Many virtualization services can redirect the USB device from the client workstation to the server/virtual desktop. PowerMic bandwidth consumption via USB redirection is approximately 440 kbit/s in Citrix environments and approximately 1 Mbit/s in other virtual environments. There can be limitations depending on the virtualization platform and version used by a site.

Note: We don't recommend using USB redirection because it depends on the state of the network and is therefore unreliable. In addition, some devices require a high network throughput even when they're idle. Audio devices work better with the native audio channel.

Use native audio channels and device splitting in the following cases:

- When using a PowerMic with a zero client or with a thin client not supported by the Nuance virtual extensions.
- When using a double-hop configuration.

Recommendations

Note: More than one Nuance solution might be needed to fulfill a given site's requirements.

Zero clients

Alternative solution: [Vendor native audio channel](#).

- Should only be considered if the site needs to use PowerMic or composite devices from other microphone vendors.
- Doesn't need additional software to be installed on the client workstation.
- The audio is compressed to optimize network bandwidth and encrypted for security (if enabled).

Remarks

- Nuance virtual extensions can't be installed on zero clients.

Single-hop configuration (thin or thick clients)

Recommended solutions:

- [Nuance virtual extensions](#) if the site uses a supported virtualization vendor.
- Vendor native audio channel if the virtualization vendor isn't supported by the Nuance virtual extensions.

Double-hop configuration (thin or thick clients)

Recommended solutions: Combination of vendor native audio channel and Nuance virtual extensions.

Recommendations overview

Technology	Bandwidth	Recommended for
Nuance virtual extensions	27 kbit/s	Microsoft Windows and Windows Embedded operating systems Azure Virtual Desktop (AVD)
Citrix HDX	Approximately 150 kbit/s	Microsoft Windows operating systems Thin and zero clients Note: We recommend setting the audio quality option to High for optimal speech recognition quality.
OmniSSA RTAV	> 200 kbit/s (depends on the OmniSSA settings)	Microsoft Windows operating systems Thin and zero clients
Citrix USB audio	500 – 1000 kbit/s	Microsoft Windows operating systems Thin and zero clients
OmniSSA USB audio	500 – 1000 kbit/s	Microsoft Windows operating systems Thin and zero clients

References

For more information, see the following:

- *Nuance - Deployment and Configuration for Citrix Environments*
- *Nuance - Deployment and Configuration for Omnissa Horizon Environments*
- *Nuance - Deployment and Configuration for Microsoft RDS Environments*
- [IGEL Knowledge Base](#)